

**DEKA**  
Innate Ability

THE NEW **PRO SYSTEM**

**NEW!**  
—  
**25 kW  
POWER**

**again**  
**PRO  
PLUS**

## **More POWER for More RESULTS**

Unmatched Performance  
for Hair Removal, Vascular  
& most Aesthetic Applications

AGAIN PRO PLUS stands out in **Aesthetics** and **Dermatology**



# AGAIN PRO PLUS

EVOLVING EXCELLENCE

HAIR REMOVAL

VASCULAR TREATMENTS

BENIGN PIGMENTED  
LESIONS

NON-ABLATIVE  
SKIN REJUVENATION

LEG VEINS

PSEUDOFOLLICULITIS  
BARBAE (PFB)

ONYCHOMYCOSIS

TOTAL FACIAL TREATMENT  
(MOVEO GLO)



# Key Principles

After the success of the AGAIN PRO platform, DEKA has designed and realized an improved platform in its style and performance: **AGAIN PRO PLUS** - with increased power levels, up to 25 kW unparalleled treatment speed and specific pulses to treat also fair and thin hair in order to maximize the efficacy of vascular treatments.

## Why Choosing **AGAIN PRO PLUS**:

### All Hair Types Are Treated


**AGAIN PRO PLUS** treats not only dark hair, but thanks to the very short pulses it is able to target also those fairer and thinner hair and in delicate areas.

### Fast, Safe and Comfortable

Treatment duration is reduced by 50% thanks to a maximized performance and an increased peak power. Thanks to a stronger cooling power, perception of pain and possible side effects are close to zero.

### Complete Multidisciplinarity in all Body Area

A flexible dual-wavelength system with a wide range of spot sizes, fluences and repetition rates; it delivers reliable, cost-effective treatments from permanent hair removal, to vascular treatments and non-ablative rejuvenation in all body areas.

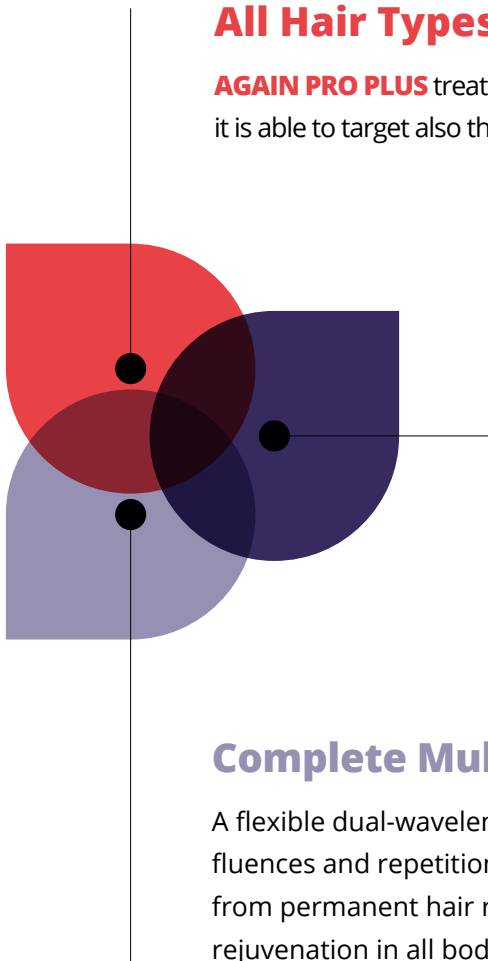



**AGAIN PRO PLUS IS SUSTAINABLE**

DEKA systems and devices are all designed and assembled to achieve our sustainability goals and to make our supply chain traceable, transparent and reliable. We match international environmental standards for ourselves and our suppliers


THE SCIENCE BEHIND

## Hair Removal Treatment: Precision, Speed and Power

**AGAIN PRO PLUS** introduces a more effective and extremely flexible dual-wavelength system delivering reliable, cost-effective treatments to all hair type (even the finest and lightest) in the fastest and safest way.


Again Pro Plus introduces Ultra Short Pulse, a cutting-edge hair removal treatment that utilizes high-intensity laser pulses to target hair follicles with exceptional precision. This advanced technology ensures fast, effective, results, with minimal discomfort.

Treatment time is reduced by 50% and performance is increased thanks to the wide range of spot sizes and higher system power.

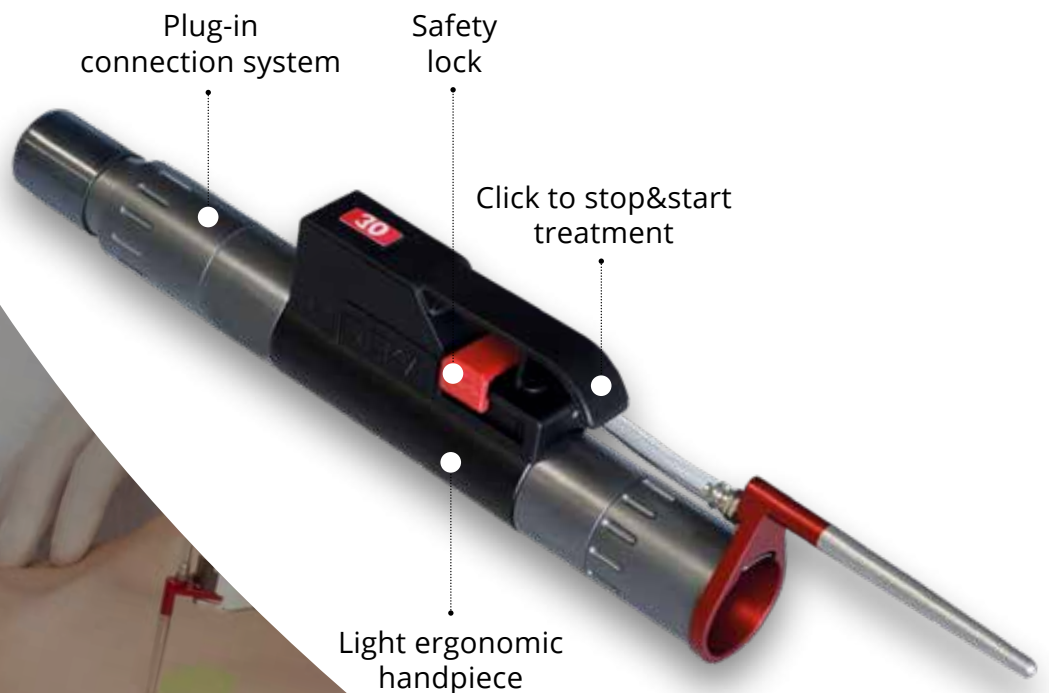



- More than 11 selectable handpieces
- Different spotsizes available (2.5-30 mm)
- Automatic spot recognition system
- Ergonomic connector for air-cooling devices
- The biggest 30 mm spotsize

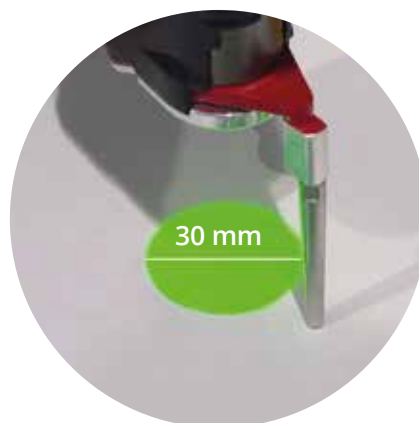

## VASCULAR TREATMENT: LONG PULSE 1064 nm

**AGAIN PRO PLUS** exploits an optimized pulse for vascular treatments which maximizes efficacy while reducing pain.



THE SCIENCE BEHIND

## Inside Moveo Innovation

**Moveo handpieces** drastically reduces the amount of energy reflected by the skin by working in contact with the skin at a lower fluence. This exclusive technology by DEKA is suitable for a wide range of applications.

moveo **HR**

**Focus on Hair Removal**

moveo **PL**

**Focus on Pigmented Lesions**

moveo **SR**

**Focus on Skin Rejuvenation**

moveo **VL**

**Focus on Vascular Lesions**



moveo**Glo**

### GLO'ing Results with Minimal Downtime

**3** Skin Face Enhancing Treatments  
in **ONE** Session

**SR** + **VL** + **PL**

With the Moveo Glo stacked treatment, AGAIN PRO can address vascular lesions, collagen stimulation, and pigmented lesions with the comfort of moveo and minimal downtime.





## Clinical Results




Courtesy of  
Prof. P. Bonan, M.D. and M. Troiano, M.D.  
*Donatello Clinic, Florence, Italy*



Courtesy of  
Prof. P. Bonan, M.D. and M. Troiano, M.D.  
*Donatello Clinic, Florence, Italy*



Courtesy of  
Setsuko Clinic,  
*Montreal, Canada*

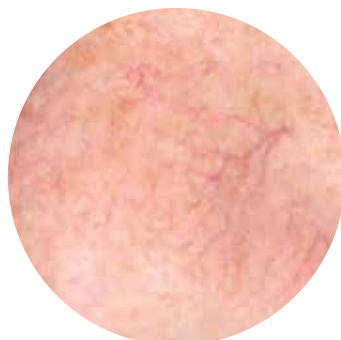

Courtesy of  
Prof. G. Cannarozzo M.D.  
*University of Rome Tor Vergata,  
Roma, Italy*

# Knowing Our Audience and Market

Thanks to growing beauty consciousness among consumers, worldwide, the global hair removal market is witnessing a high-paced increase in its size and valuation.

## Global Hair Removal Market\*



**1st**  
**REQUESTED TREATMENT**  
within the Aesthetic Industry



**+17%**  
**COMPOUND ANNUAL  
GROWTH RATE**



**↑ 3.9** Billion \$  
**EXPECTED UNTIL 2026**

**\*Sources:**

Market Study Report "Global Laser Hair Removal Market Size Status And Forecast 2019-2025"

Acumen Research and Consulting "Laser Hair Removal Market Size, Share, Trends, Scope, Growth Opportunity and Forecast 2019-2026"





## The Floor to Practitioners

“ AGAIN PRO PLUS is the latest generation laser device and it's the evolution of AGAIN PRO but it's recently been improved in the design, user interface and it's been further enhanced with a more powerful and faster Alexandrite together with a more powerful and, at the same time, very gentle Nd:YAG.

For vascular application the Nd:YAG 1064 has now the capability to be extremely delicate, to be extremely gentle thanks to the possibility to split the pulses: The first one which is able to convert the hemoglobin into metahemoglobin allowing the second one to absorb this target more thereby closing the vessel.

The cooling power of Moveo for AGAIN PRO PLUS has improved by 40% allowing for the treatment to be more gentle and safer for darker phototypes and has been improved with an increased pulse emission frequency which allows for the treatment of patients in a very fast manner. The 30mm spotsize allows the operator to perform hair removal treatments quickly, so we can complete a total-body procedure in less than 20 minutes and the back of a male in less than 2 minutes by using 2.5 or 3 Hz, thanks to the increased power of the device.

Moreover, the 30mm spotsize is able to penetrate very deep into the tissue, about 5 mm and this way it can even treat very deep masculine hairs, an ever-growing target for us.

**Dr. Michela Troiano, MD- Dermatologist**  
*(Italy)*

THE PERFECT MATCH

## AGAIN PRO Strength

- **Up to 25 kW power** for hair Removal treatments.
- **High-powered laser technology**  
flexibility with a wide range of spot sizes,  
fluences and repetition rates for fast and efficient treatments.
- **A modern 2-in-1 workstation**  
performing the most requested aesthetic applications  
(hair removal, vascular lesions and skin spots, non-ablative  
photorejuvenation).
- Suitable **for all skin types**.
- **Fast, safe & painless** with Moveo handpieces.
- **Easy switch between wavelengths**  
and different treatments within a single session.
- **Know-how at your fingertips**  
with the built-in protocol database for all kinds of hair,  
phototype & treatment.
- **Low maintenance cost.**



**Only for DEKA Users:**

A great support in medical practice

**DEKA** Club

**Extended tube for external cooling system**

**180° rotating and tilting LCD 15.6"**

**Vascular pulses for Nd:YAG**

**Optical fiber**

**More than 11 standard handpieces (2.5-30 mm)**

**Optional Moveo technology (HR, PL, SR, VL)**

**Alex & Nd:YAG Double source**



# Technical Specifications

AGAIN PRO PLUS		
Laser Type	755 nm Alexandrite*	1064 Nd:YAG
Energy per Pulse	64 J (max)	120 J (max)
Laser Power	125 W (max)	180 W (max)
Pulse Duration	0.2 – 300 ms	
Repetition Rate	Up to 12 Hz	
Spot Sizes (optional)	From 2.5 mm to 30 mm, Moveo HR, Moveo PL, Moveo SR, Moveo VL	
Fluence	Up to 600 J/cm <sup>2</sup>	
Dimensions and Weight	130 (H) cm x 40 (W) cm x 105 (D) cm   150 Kg	
Electrical Requirements	<30A @230 Vac	

\* A single wavelength version is available, only equipped with Alexandrite laser

**DANGER** - Visible and invisible laser radiation. Avoid eye or skin exposure to direct or scattered radiation. Class 4 laser product.

This brochure is not intended for the US market.



Follow us on



[www.dekalaser.com](http://www.dekalaser.com)

Dealer stamp



**DEKA M.E.L.A. s.r.l.**

Via Baldanzese, 17 - 50041 Calenzano (FI) - Italy  
Tel. +39 055 8874942 - Fax +39 055 8832884

#### DEKA Innate Ability

A spin-off of the EL.En. Group, DEKA is a world-class leader in the design and manufacture of lasers and light sources for applications in the medical field. DEKA markets its devices in more than 120 countries throughout an extensive network of international distributors as well as direct offices in Italy, France, Japan and the USA. Excellence is the hallmark of DEKA's experience and recognition garnered in the sphere of R&D in over thirty years of activity. Quality, innovation and technological excellence place DEKA and its products in a unique and distinguished position in the global arena. DEKA manufactures laser devices in compliance with the specifications of Directive 93/42/EEC and its quality assurance system is in accordance with the ISO 9001 and ISO 13485 standards.\*